

# danumed® Extension Set for Gastrostomy Button

**Security  
meets Comfort**

**Durable TRITAN™  
Connector Material ✓**

**Kink-resistant  
Polyurethane Tube ✓**

**Innovative Rotating  
Ring ENSwivel® ✓**



**German  
Manufacturer ✓**

**Reliable  
Supply Chain ✓**

**Optimised Button Connectors,  
Straight and Angled ✓**


**Long Application Time  
up to 14 Days ✓**



# danumed® Extension Set

The danumed® Extension Set is designed for the administration of enteral nutrition, liquid medication and for decompression in connection with a Gastrostomy Button.

Innovative features optimise the flow rate and simplify the handling for the user.

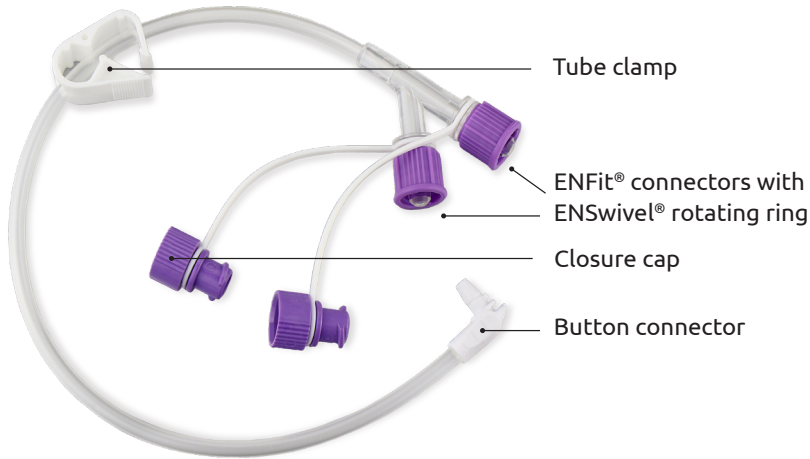

**High-quality materials allow for up to 14 days of use:**

## TRITAN™ Connector Material


- extremely resistant to aggressive drugs and solutions
- break-resistant, can avoid stress cracks

## Polyurethane Tube Material

- kink-resistant
- reduces the hazard of tube occlusions
- biocompatible, without harmful substances


## Innovative Rotating Ring ENSwivel®

- For enhanced patient comfort and easy user handling
- Easy connection and release of even very tight ENFit® connections
- Twisting of tubes is avoided
- Fitted on all ENFit® connectors

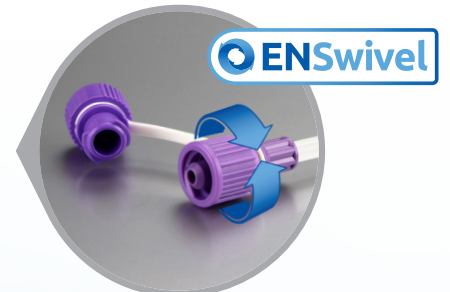


## High-End Button Connector

- Rounded inner profile optimises the flow rate, independent of the tube lumen
- Easy, secure and tight connection to the Button

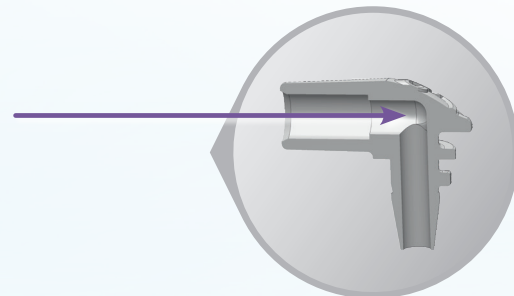

## Numerous Versions and Lengths

- Available ENFit® connector versions: single ENFit®, Y-ENFit®, Bolus-ENFit®
- Available in different lengths: 5 cm, 30 cm and 60 cm

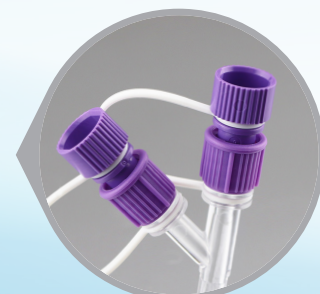




### danumed® Extension Set with angled Button Connector

REF	Description	Length	Feeding connector
DA 1011 3011	Extension Set <b>sterile</b> , 30 cm, ENFit®, ENSwivel®	30 cm	ENFit®
DA 1010 3011	Extension Set <b>non-sterile</b> , 30 cm, ENFit®, ENSwivel®	30 cm	ENFit®
DA 1011 6011	Extension Set <b>sterile</b> , 60 cm, ENFit®, ENSwivel®	60 cm	ENFit®
DA 1010 6011	Extension Set <b>non-sterile</b> , 60 cm, ENFit®, ENSwivel®	60 cm	ENFit®
DA 1011 0521	Extension Set <b>sterile</b> , 5 cm, Y-ENFit®, ENSwivel®	5 cm	Y-ENFit®
DA 1011 3021	Extension Set <b>sterile</b> , 30cm, Y-ENFit®, ENSwivel®	30 cm	Y-ENFit®
DA 1010 3021	Extension Set <b>non-sterile</b> , 30 cm, Y-ENFit®, ENSwivel®	30 cm	Y-ENFit®
DA 1011 6021	Extension Set <b>sterile</b> , 60 cm, Y-ENFit®, ENSwivel®	60 cm	Y-ENFit®
DA 1010 6021	Extension Set <b>non-sterile</b> , 60cm, Y-ENFit®, ENSwivel®	60 cm	Y-ENFit®



Delivered with right-angled button connector and **standard bore tube**. Individually packaged. 10 pieces in a box, 10 boxes in a shipping box.

## Product Versions

### Bolus Nutrition with Syringe

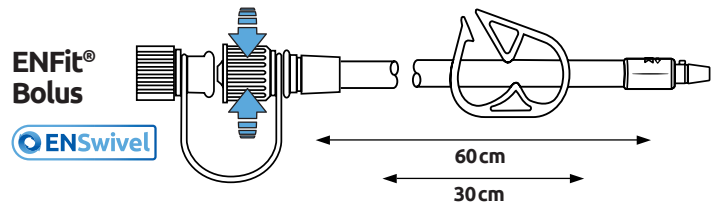

### Continuous Nutrition with Pump or Gravity Set

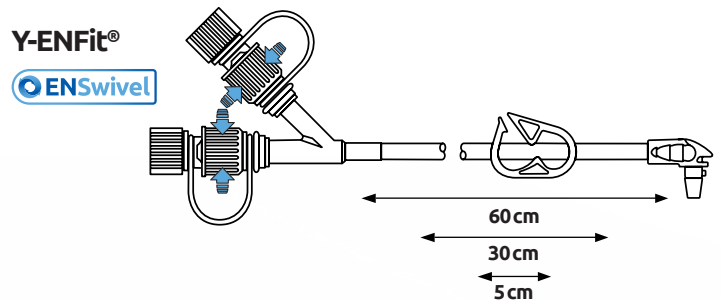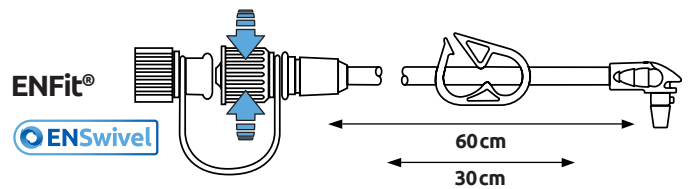

### danumed® Extension Set with straight Button Connector

REF	Description	Length	Feeding connector
DA 1011 3032	Extension Set <b>sterile</b> , 30 cm, Bolus-ENFit®, ENSwivel®	30 cm	Bolus-ENFit®
DA 1010 3032	Extension Set <b>non-sterile</b> , 30 cm, Bolus-ENFit®, ENSwivel®	30 cm	Bolus-ENFit®
DA 1011 6032	Extension Set <b>sterile</b> , 60 cm, Bolus-ENFit®, ENSwivel®	60 cm	Bolus-ENFit®
DA 1010 6032	Extension Set <b>non-sterile</b> , 60 cm, Bolus-ENFit®, ENSwivel®	60 cm	Bolus-ENFit®

Delivered with straight button connector and **wide bore tube**. Individually packaged. 10 pieces in a box, 10 boxes in a shipping box. The Bolus Extension Sets are available on request with funnel connector instead of ENFit® connector.



## **danumed Medizintechnik GmbH**

*Individual solutions are our standard.*

We are a family-run, innovative, German medical technology company with high quality standards, many years of expertise in enteral nutrition and a clear customer focus.

Individual customer requests are met with flexibility and reliability.  
Personal, competent, cooperative.

For more information, please contact us on **(02) 9882 3666** or at **[amsl@amsl.com.au](mailto:amsl@amsl.com.au)**

**[amsl.com.au](http://amsl.com.au)**

ALWAYS READ THE LABEL AND FOLLOW THE DIRECTIONS FOR USE.  
Manufactured by danumed. Distributed in Australia by AMSL.  
AMSL is a Dexcom company. PR-100-764 July 2023

