

SE-601 Series

Electrocardiograph



SE-601 Series

Electrocardiograph



About Edan

Edan is a healthcare company dedicated to improving the human condition around the world by delivering value-driven, innovative and high-quality medical products and services. For over 20 years, Edan has been pioneering a comprehensive line of medical solutions that address a broad range of healthcare practices including:

- Diagnostic ECG
- Patient Monitoring
- OB/GYN
- Ultrasound Imaging
- Point-of-Care Testing
- *In-Vitro* Diagnostics
- Veterinary

Healthcare professionals around the world depend on Edan's breakthrough medical technologies and outstanding customer support.



Global Headquarters:
Edan Instruments, Inc. | 15 Jinhui Road, Pingshan District, Shenzhen
518122 P.R. China | +86.755.26898326 | www.edan.com | info@edan.com

U.S. and Canada inquiries:
EDAN Diagnostics, Inc. | 9918 Via Pasar, San Diego, CA 92126
+1.858.750.3066 | www.edandiagnostics.com | edan-info@edandiagnostics.com

© Edan Instruments, Inc. All rights reserved. Features and specifications are subject to change without prior notice. No reproduction, copy or transmission may be made without written permission. Not all products or features are available in all countries, contact Edan for local availability.



SE-601 Series

Multi-channel Electrocardiograph



User-friendly Design

- Compact and portable design
- 5.6" high resolution color touch screen(optional)
- Alphanumeric keyboard

Streamlined Workflow

- One-touch operation
- Support barcode scanner to simplify workflow (for SE-601B/C)

Accurate Clinical Performance


- 12 leads simultaneous acquisition
- Sampling rate as high as 16,000 Hz aims at capturing all the details to generate precise result
- Weak signal detection with the advanced 24-bit A/D converter technology
- Anti-noise technology with comprehensive filters and the high CMRR over 140dB
- Bandwidth from 0.01 to 300 Hz guarantees the interpretation especially for pediatric patients
- Extend printing once arrhythmia is detected to further your rhythm analysis
- Automatic measurements and interpretation tested with CSE/AHA/MIT database
- Optional Glasgow Algorithm with specialized Children ECG Interpretation

Data Management

- Internal storage up to 800 ECGs and can be enlarged by USB flash disk
- SE-1515 software for data management on PC (optional)
- Report formats: PDF(for SE-601B/C) and optional
- SCP/FDA-XML/DICOM data export



SE-601C With Rolling Stand
(5.6" Color Touch Screen)

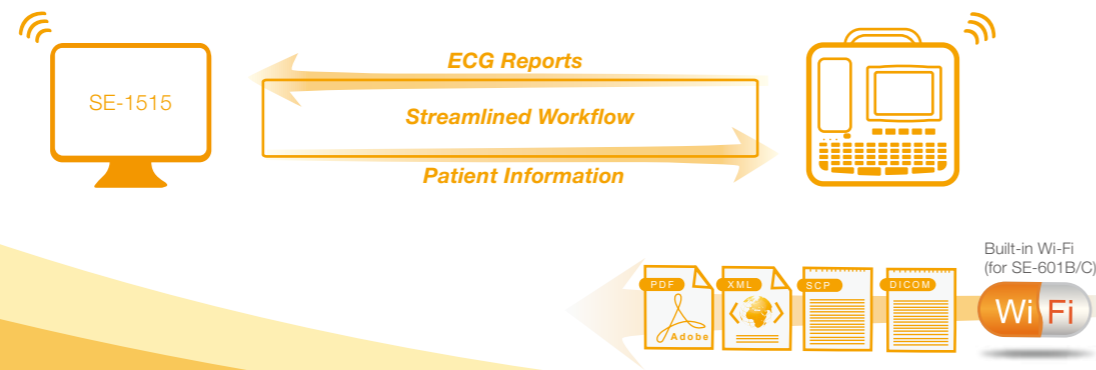

SE-601A
(3.5" Monochrome LCD Display)



SE-601B
(5.6" Color LCD Display)

Seamless Connectivity



Network your ECG data with the help of EDAN SE-1515 Data Management System

- Bi-directional communication with ECG machines on worklist and reports
- Cardiac consultation upon ECG reports over database sharing
- Multiple report formats output including PDF, XML, SCP, DICOM, JPG, BMP, and DOC
- Data access with HIS/EMR/PACs/CIS based on HL7/DICOM/GDT protocol