

Suprasorb® P sensitive

Versatile, protective
and gentle.



Suprasorb® P sensitive

Five layers of softness for patient comfort.

Can be used as part of
a prophylactic therapy
to support the prevention
of pressure ulcers.*

Protective backing:

- High moisture vapour transmission rate (MVTR)
- Smooth and conformable
- Waterproof
- Bacterial barrier

Suprasorb P sensitive border with super-absorbent core:

- Super-absorbent polymer fibres
- Good absorption capacity
- Helps reduce the risk of maceration

Non-woven layer:


- Helps to distribute exudate

Absorbent PU foam*:

- Efficient exudate management
- Offers comfort and a cushioning effect


OptiSil silicone wound contact layer:

- Technology combining perforations and silicone oil
- Enables quick absorption of exudate into the dressing
- Soft and minimises pain during dressing removal
- Limits the risk of adhesion to the wound

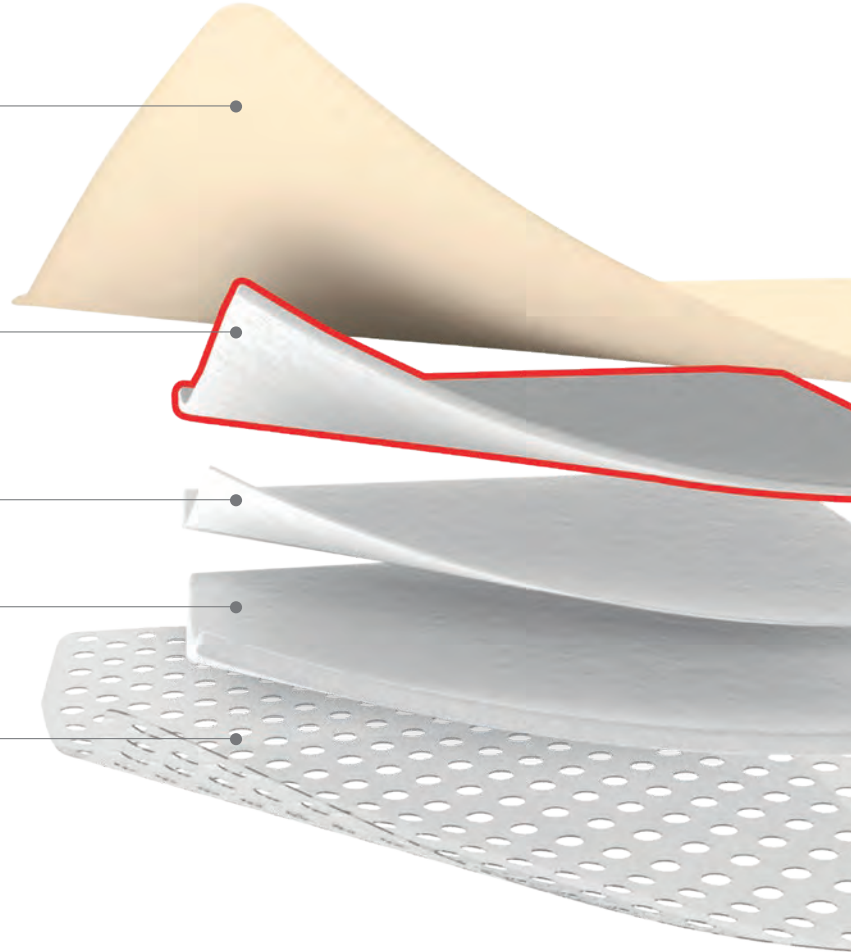

* Polyurethane foam

Easy application:

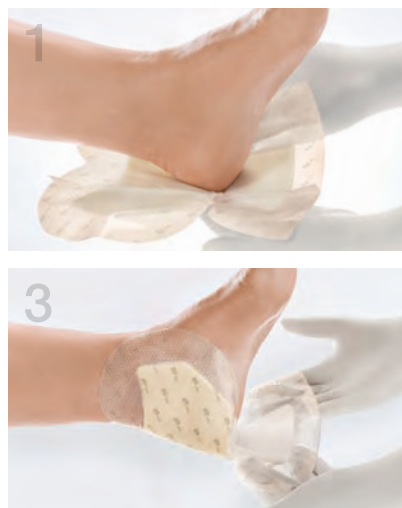

With 2 or 3 liners, which helps to avoid touching the adhesive with gloves.

Suprasorb P sensitive is versatile and can even be used for wounds in **awkward areas**.

Suprasorb® P sensitive

Versatile wound care even with fragile skin.



Protection

- water- and bacteria-proof
- permeable to water vapour

Efficient exudate management

- promotes quick absorption of exudate into the dressing
- maintains moist environment



Softness & comfort

- adheres gently
- limits the risk of adhesion to the wound during dressing changes
- comfortable thanks to absorbent PU foam



Silicone wound contact layer suitable for **fragile skin**



Efficient exudate management reduces the **risk of maceration** and promotes wound healing



Can be worn **in the shower**



Minimises **pain** to the patient and **trauma** to the wound and skin during dressing removal or repositioning

The new Suprasorb P sensitive is a next-generation silicone foam dressing for a broad range of wound types. This innovative wound dressing provides numerous benefits for the patient. Up to five different material layers provide efficient exudate management to promote wound healing while protecting the skin.

Available in 3 versions:

border - 5 layers

with super-absorbent core




border lite - 4 layers

without super-absorbent core



non-border - 3 layers

without super-absorbent core

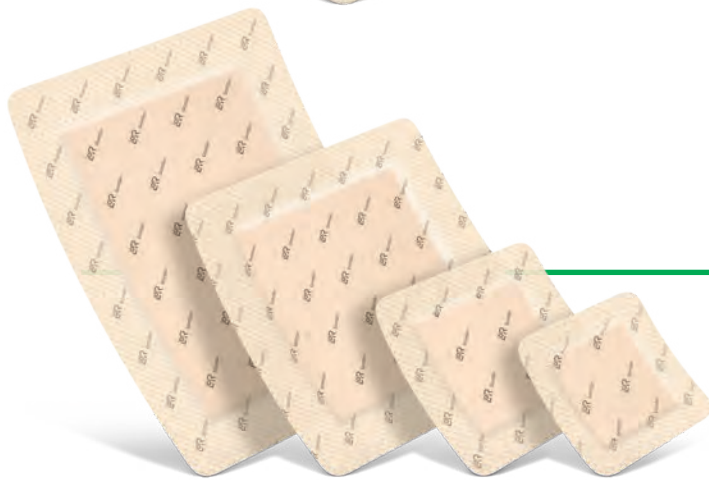

The right shape for every task.

The products **border**, **border lite** and **non-border** are available in different formats and sizes.



Suprasorb® P sensitive tracheo - 3 layers

without super-absorbent core
For tracheostoma wound care.



Suprasorb® P sensitive sacrum - 5 layers

with super-absorbent core
Anatomically shaped for use around the sacrum area.



Suprasorb® P sensitive multisite - 4 & 5 layers

with and without super-absorbent core
For the care of awkward wound areas.



Suprasorb® P sensitive heel - 5 layers

with super-absorbent core
For wound care on joints and heels.

Suprasorb[®] P sensitive

Dressings for a broad spectrum of wounds.

border
5 layers



for moderately to heavily exuding wounds

border lite
4 layers



for lightly to moderately exuding wounds

non-border
3 layers




for moderately exuding wounds



Superficial wounds



Acute wounds

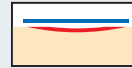


Chronic wounds



Deep wounds*

* In combination with a suitable primary dressing.

Indicated for superficial to deep wounds (in combination with a suitable primary dressing) such as:

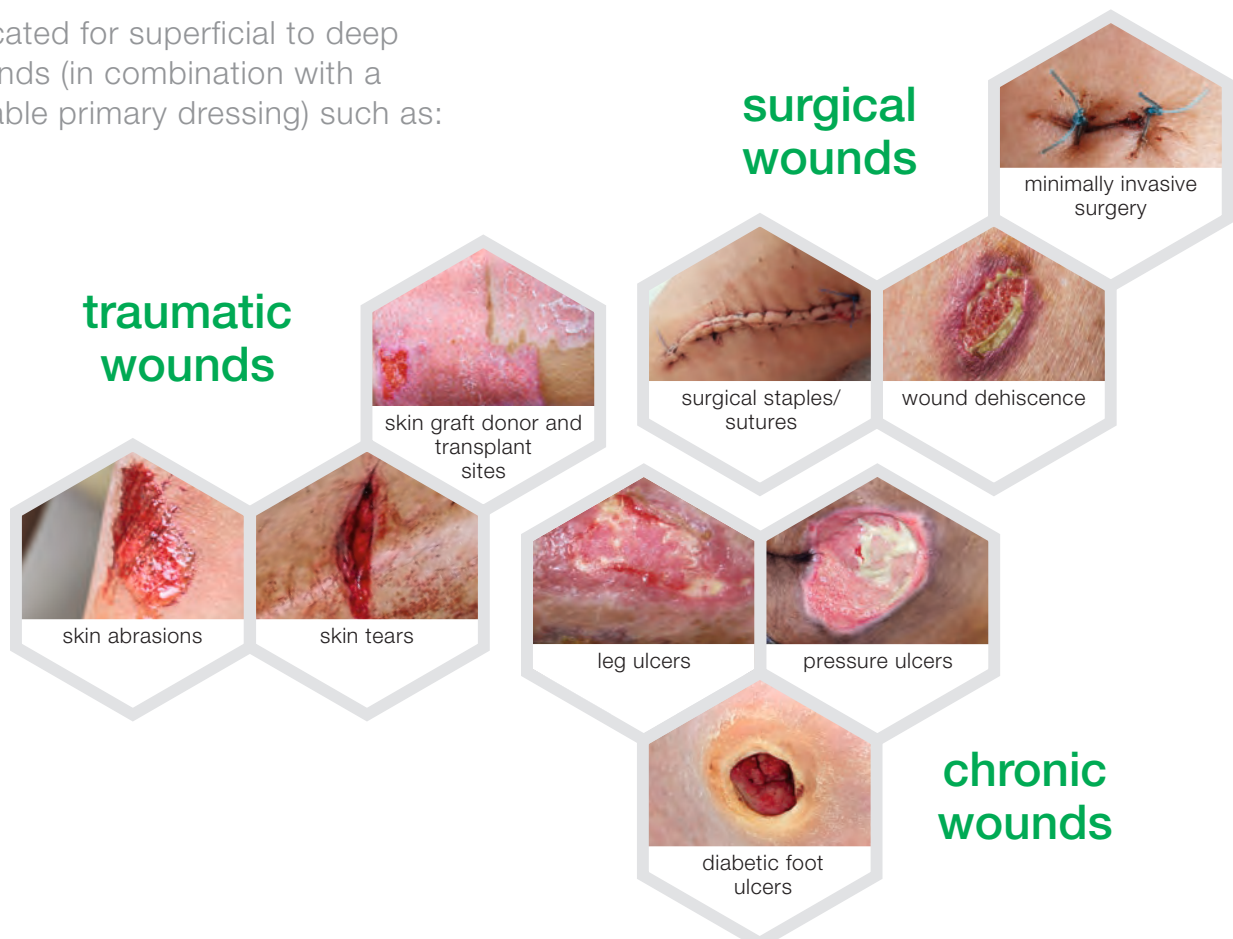




Efficient exudate management

Protection

Softness & comfort

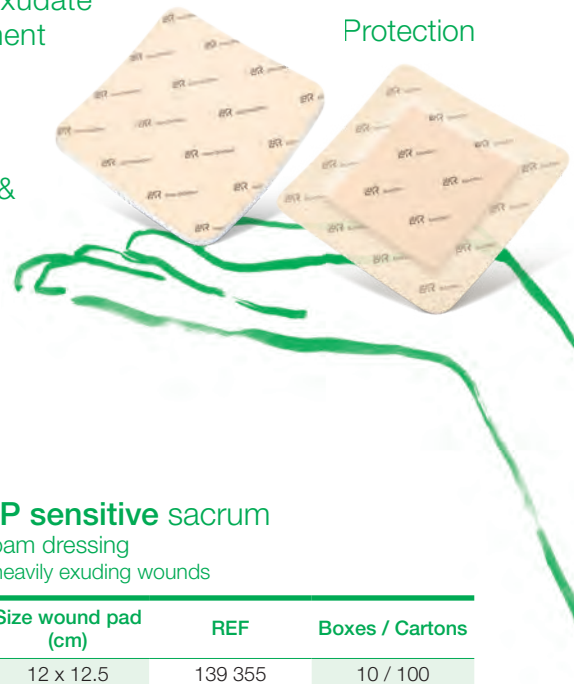

Suprasorb® P sensitive border

5-layer silicone foam dressing
For moderately to heavily exuding wounds

Size dressing (cm)	Size wound pad (cm)	REF	Boxes / Cartons
7.5 x 8.5	4.5 x 5	139 068	10 / 200
10 x 10	6.2 x 6.2	139 061	10 / 200
12.5 x 12.5	8.5 x 8.5	168 062	10 / 200
15 x 15	10.5 x 10.5	168 069	10 / 150
10 x 30	24 x 5	146 240	5 / 200
20 x 20	15 x 15	139 070	10 / 100

Suprasorb® P sensitive border lite

4-layer silicone foam dressing
For lightly to moderately exuding wounds

Size dressing (cm)	Size wound pad (cm)	REF	Boxes / Cartons
5 x 5	2.5 x 2.5	139 080	10 / 240
5 x 10	2.5 x 6	139 081	10 / 200
7.5 x 8.5	4.5 x 5	139 086	10 / 200
9 x 15	5 x 10.5	168 087	10 / 200
10 x 10	6.2 x 6.2	139 083	10 / 200
10 x 20	5 x 15	139 088	10 / 150

Suprasorb® P sensitive tracheo

3-layer silicone foam dressing
For moderately exuding wounds

Size dressing (cm)	Size wound pad (cm)	REF	Boxes / Cartons
10 x 9	9 x 10	146 244	10 / 100

Suprasorb® P sensitive non-border

3-layer silicone foam dressing
For moderately exuding wounds

Size dressing (cm)	REF	Boxes / Cartons
5 x 5	139 350	10 / 240
7.5 x 7.5	139 455	10 / 150
10 x 10	139 351	10 / 100
10 x 20	139 352	10 / 150
15 x 15	139 353	10 / 100
15 x 20	139 456	10 / 100
20 x 20	139 354	10 / 100

Suprasorb® P sensitive sacrum

5-layer silicone foam dressing
For moderately to heavily exuding wounds

Size dressing (cm)	Size wound pad (cm)	REF	Boxes / Cartons
17 x 17.5	12 x 12.5	139 355	10 / 100
23 x 23	16.5 x 17	139 556	10 / 100

Suprasorb® P sensitive heel

5-layer silicone foam dressing
For moderately to heavily exuding wounds

Size dressing (cm)	Size wound pad (cm)	REF	Boxes / Cartons
23.5 x 25	17.5 x 19	139 356	10 / 100

Suprasorb® P sensitive multisite border

5-layer silicone foam dressing
For moderately to heavily exuding wounds

Modified design, high conformability

Size dressing (cm)	Size wound pad (cm)	REF	Boxes / Cartons
7.5 x 9.5	3.5 x 5.5	168 075	10 / 300
12 x 15	7 x 10	168 357	10 / 100
16 x 20	11.5 x 15.5	168 076	10 / 200

Suprasorb® P sensitive multisite border lite

4-layer silicone foam dressing
For lightly to moderately exuding wounds

Modified design, high conformability

Size dressing (cm)	Size wound pad (cm)	REF	Boxes / Cartons
7.5 x 9.5	3.5 x 5.5	168 077	10 / 300

For any queries please contact your local Lohmann & Rauscher (L&R) Sales Representative or L&R Customer Service on:

Australia: 1300 572 869 New Zealand: 0800 572 869 Email: info@au.LRmed.com

* World Union Of Wound Healing Societies Consensus Document (WUWHS) - Role of dressings in pressure ulcer prevention, Box 5 | Ideal properties of a dressing used for PU prevention page 11. L&R Data on file demonstrates alignment with relevant properties as referenced in the Consensus Document.

Suprasorb® is a registered trademark of Lohmann & Rauscher GmbH & Co KG. © 2024